

# PAUSE M4

---

USER MANUAL



**CERAGEM**



Thank you for purchasing the CERAGEM Massage Chair.  
Please read this manual carefully before starting any use.  
Please store the manual where you can always find it after reading it.

This User's Manual contains the product warranty information.  
This product is a home massage chair for household use.  
Do not use for commercial or industrial purposes.  
We may improve the design and specifications of this product  
without any further notice.

This product is exclusively for 110V.  
Please make sure to plug the power plug into a dedicated 110V outlet.  
Do NOT use a 220V transceiver to 110V as it may result in safety hazards.



Users whose height exceed 6'3" or weight more than 220 lbs is  
not recommended to use the product. Product may fall over and  
get damaged which may cause accidents or injuries to user.


# CERAGEM LIFESTYLE MASSAGE CHAIR USER MANUAL

PAUSE M4 USER MANUAL

MASSAGE CHAIR  
CGMEMCG-2301

# CONTENTS

- 01 Safety Precautions
- 02 Care and Storage
- 03 Assembly Method
- 04 Name of each part and the components
- 05 How to use the remote controller
- 06 Troubleshooting
- 07 FAQ
- 08 Customer Service
- 09 Product Specification
- 10 Product Warranty
- 11 FCC Warning

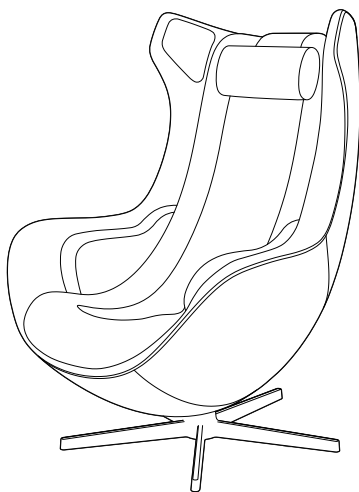

**MASSAGE CHAIR**  
**CGMEMCG-2301**

## IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this appliance:

**DANGER** – To reduce the risk of electric shock:

- 1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.

**WARNING** – To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- 4) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 5) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- 6) Do not carry this appliance by supply cord or use cord as a handle.
- 7) Keep the cord away from heated surfaces.
- 8) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 9) Never drop or insert any object into any opening.
- 10) Do not use outdoors.
- 11) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- 12) To disconnect, turn all controls to the off position, then remove plug from outlet.
- 13) Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- 14) Connect this appliance to a properly grounded outlet only. See Grounding Instructions.

The following safety precautions are intended to prevent unexpected danger or damage by ensuring safe and accurate use of the product.

### **When Installing**

Remodeling, disassembly, or repairs should not be performed by anyone other than service personnel.

→ Can cause fire or electric shock.

Use the product on a flat surface.

→ Can cause product failure or an accident.

### **Regarding Power**

If the product does not work or ceases to function, immediately disconnect the power plug from the outlet and contact the customer center for inspection and repair.

→ There is a risk of fire or electrocution due to electric shock or short circuiting.

If the cover of the product is torn or peeled off, it may cause a serious accident as clothes or hair can be caught inside. Ensure that there are no issues with the appearance of the product. If there is any damage, please contact the support center.

→ It may cause electric shock or injury.

Instruct children not to touch the power plug.

→ There is a risk of electric shock or injury.

Do not bend the power cord too hard and do not place heavy objects on it.

→ There is a risk of fire and electric shock.

Do not use multiple devices simultaneously when using a multi outlet adapter.

→ There is a risk of fire due to overheating of the outlet.

When unplugging the plug, pull it out of the outlet by grasping the plug.

→ If the power cord is severed, there is a risk of fire and electric shock.

Do not touch the power plug with wet hands, and when you clean or repair the product, disconnect the power plug from the outlet.

→ There is a risk of electric shock.

This product is exclusively for 110V. Please make sure to plug the power plug into a dedicated 110V outlet.

→ There is a risk of fire due to overheating of the outlet

Be sure to use a grounding outlet when using an extension cord.

→ Using a non grounding outlet may cause electric shock or injury.

When there is a power outage, immediately turn off the power switch and unplug the power plug from the outlet.

→ Can cause electric shock or fire.

While operating the product, avoid pinching the AC cord or remote controller line between the center connection support or the five foot support.

→ A short circuit may occur due to damage to the covering, so be careful when using it.

### **When using the product**

Due to the risk of safety accidents, users with the following conditions should refrain from using the product:

The following persons should not use this product:

- People with thrombosis, embolism, severe venous, acute venous, etc.
- People with suspected acute cervical sprain (whiplash)
- People with symptoms of a sprain or sudden muscle contraction
- Pregnant women in early pregnancy or immediately after childbirth
- Those who have an abnormality in their spine due to an accident or disease in the past, or those with bent spine (scoliosis)
- Those who have an acute disease of the intestines or severe heart disorder
- Those who use an implantable electronic medical device such as implantable pacemakers
- Those who have a severe seizure symptom or psychogenic disorder (perceptual disorder)

- 
- Those whose temperature sensation is lost due to severe peripheral circulation disorders such as diabetes, etc.
  - Those who have a risk of osteomyelitis or rickets
  - \*Rickets: Symptoms of shrinking muscles throughout the body and an inability to stretch arms and legs
  - Those who are unable move their body freely
  - Those with or undergoing treatment for symptoms of disc disorders
  - Those receiving treatment at a medical institution for any disorder in their body

The following persons must consult a doctor before use:

- Those who have a risk of side effects due to the use of the product, e.g. pregnant women
- Those whose body temperature is higher than normal due to acute inflammation, chills, etc.
- Those with symptoms of acute skin disease, malignant tumor, or other skin disease
- Those who are receiving medical treatment
- Those who feel any bodily abnormalities
- Those who have been instructed by doctor not to undergo massage
- Those whose muscles have weakened or become lean due to old age
- Those who suffer from back pain

Do not use or store the product in a humid place such as a bathroom.

- It may cause electric shock, fire, failure, or mold.

Do not allow children or pets to play near or on top of the product and prevent them from using it.

- It may cause malfunction or injury.

Use under supervision if the user has a cognitive impairment or is inexperienced with the product.

- It may cause malfunction or injury.

Do not wear something hard on your head such as hair accessories during use.

- It may cause injury.

The operation time should not exceed a total of 30 minutes per day and 15 minutes at one time. Continuous use targeting the same area should be 5 minutes or less.

- Prolonged use may strain your muscles or nerves.



If you experience any of the following side effects while using the product, stop using it:

- Abrasions, strains, bruises, stinging/burning, or cuts caused by the product
- Damage to muscles, tendons, nerves, spine, etc. due to excessive use and pressure from the product
- Skin diseases and burns such as fractures, dislocations, sprains, etc.
- Respiratory and digestive disorders including asphyxiation, etc.

Do not use the product if you have any objects inside your pockets when massaging your hips or thighs.

- There is a risk of an accident or injury.

Make sure to use the product after checking that the power and remote control cords are neatly organized.

- It may cause an accident or injury or product failure.

Do not disconnect the power plug from the wall outlet or turn the power switch to [OFF] during the massage operation.

- This may result in an injury.

Do not use near open flames such as stoves or while smoking. Also, do not use it while mounted on or placed on top of a heating apparatus.

- It may cause fire.

Do not use after drinking alcohol.

- This may result in an accident or injury.

When operating product, please use the product by one person per product at a time. (Do not allow more than 2 or more persons to sit on the product)

- It may cause malfunction or injury.

Please do not use it on your bare skin.

- Direct massage on bare skin may cause skin damage.

Be careful not to get your legs caught in the remote control cable.

- This may result in an injury.

When cleaning or servicing the product, disconnect the power plug from the outlet.

- It may cause an electric shock or injury.

---

Do not sit while reclining.

- It may break down, or the device may fall and hurt itself.

Before the start of the massage, or through reclining, make sure there are no human (especially children) or pets around when laying down or pulling up the massage chair.

- There is a risk of accident or injury.

When using the massage chair lying down, be careful of accidents such as abalone when a person around you presses it.

- It causes an accident, injury, or product failure.

Do not use it continuously because there may be a risk of low-temperature burns on the skin if used for a long time.

- Depending on the user's condition, you can select 95 °F/104 °F/113 °F from the 3 stages of thermal intensity.

Don't let the bare skin touch the area where the heat is activated

- There is a risk of burns.
- Children with weak skin, the elderly, and patients need the protection of their guardians, and please use them more carefully.

### **When using an abdominal vibration**

When using abdominal vibration, follow the instructions and follow the usage period and method. There is a risk of burns due to incorrect usage, so please observe the following when using it.

- Use the temperature at an appropriate temperature as there may be a risk of burns if set to high temperature.
- For children, those who cannot set the temperature themselves, or those with a weak sense of skin, please be careful when using it for a long time as there is a risk of low-temperature burns despite the low temperature setting. (Make sure to wear clothes.)
- Do not use more than 2 to 3 minutes on one part of the body when setting more than 113 °F. (Please Use it for 2-3 minutes and move it to another part.)
- The surface of the product is hot, so be careful when using the product if you are sensitive to heat.
- When using a product, make sure to use a cover to prevent the product from directly touching the skin. (There's a risk of burns)


## GROUNDING INSTRUCTION

---

This product must be grounded. If it should malfunction or breakdown, grounding provides a path of least resistance for electric current to reduce the risk of electric shock. This product is equipped with a cord having an equipment grounding conductor and a grounding plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances.

**DANGER** – Improper connection of the equipment-grounding conductor can result in a risk of electric shock. Check with a qualified electrician or serviceman if you are in doubt as to whether the product is properly grounded. Do not modify the plug provided with the product – if it will not fit the outlet, have a proper outlet installed by a qualified electrician.

This product is for use on a nominal 120-volt circuit, and has a grounding plug that looks like the plug illustrated in sketch A in following figure. A temporary adapter that looks like the adapter illustrated in sketch B and C may be used to connect this plug to a 2-pole receptacle as shown in sketch B if a properly grounded outlet is not available. The temporary adapter should be used only until a properly grounded outlet (sketch A) can be installed by a qualified electrician. The green colored rigid ear, lug, or the like extending from the adapter must be connected to a permanent ground such as a properly grounded outlet box cover. Whenever the adapter is used, it must be held in place by a metal screw.

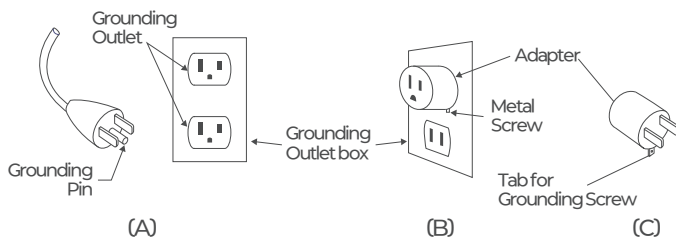

No adapter should be used with this product in Canada

---

This product must be connected to a grounded metal, permanent wiring system; or an equipment-grounding conductor must be run with the circuit conductors and connected to the equipment-grounding terminal or lead on the product.

The appliance has no user serviceable parts.

**WARNING** – This appliance is equipped with a thermal resetting protection device. This is a safety feature to protect the appliance from overheating. If the appliance should suddenly stop and will not start, turn off the main switch and pull the plug off the socket outlet. Do not operate the appliance for at least 30 minutes. Failure to turn the appliance off may result in the appliance starting unexpectedly when the device cools.

## **SAVE THESE INSTRUCTIONS**

### **DANGER**

- 1) Never use pins or other metallic fasteners with this appliance.
- 2) Carefully examine the covering before each use. Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering, or cracking.
- 3) Keep Dry - Do not operate in a wet or moist condition.

### **WARNING**

- 1) Temperature sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant or invalid or on a sleeping or unconscious person. Do not use on insensitive skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.
- 2) Do not crush - Avoid sharp folds.
- 3) Do not stand on or in the appliance. Use only while seated.
- 4) Do not use massager in close proximity to loose clothing or jewelry.
- 5) Keep long hair away from massager while in use.

When cleaning or maintaining the massage chair, be sure to disconnect the power plug from the outlet. Do not touch it with wet hands.

- There is a risk of an electric shock or injury.

Do not wipe with benzene thinner (diluent), alcohol or spray pesticides.

- It may cause discoloration, deterioration or fire.

Don't douse the product with water. Please wipe gently with a leather cleaner which is used for natural leather products.

Never iron the leather finish.

The leather finish is prone to color bleeding, so make sure to avoid wearing clothes susceptible to discoloration.

Continuous exposure to sunlight and ultraviolet rays may cause product discoloration. Avoid installation near windows or in areas with strong sunlight.

### **Information on natural leather**

PAUSE uses natural cow leather.

It is a natural phenomenon that each natural leather may vary in texture and color. Marks such as scars and wrinkles are rather characteristics of genuine leather.


Due to the features of the natural cowhide, the density of wrinkles and the texture of the skin may vary, so the surface pattern is not uniform, having slight variations.

- Leather is subject to different degrees of dyeing depending on its surface condition. The resulting slight color variation is natural appearance.


## 03 Assembly Method

---

1. Safely remove all the parts enclosed in the box.  
(Please make sure that there are no parts left behind.)
2. Place it on a flat surface to avoid scratching.
3. Lay the product forward, and then insert the center connecting support into the groove on the bottom of the product.  
\* Caution: Fit the first side of the 'center connecting support' to the fixing device.
4. Slide the "Five-leg supporting board" all the way into the coupling groove.
5. Plug the connector to the sockets and turn the fastener to the right to fix it firmly.

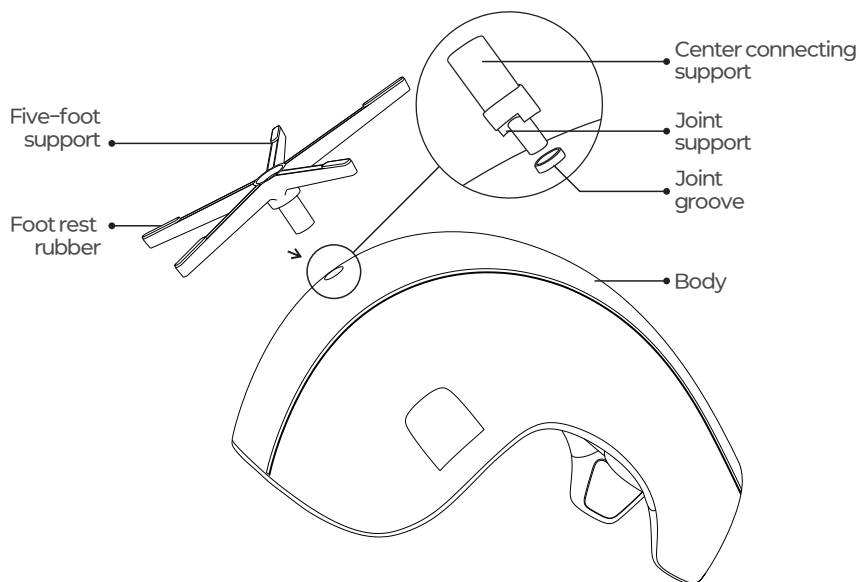


\*When installation, you must slowly drop the product.  
If you put it down suddenly, there is a risk of damage  
or injury to the product.

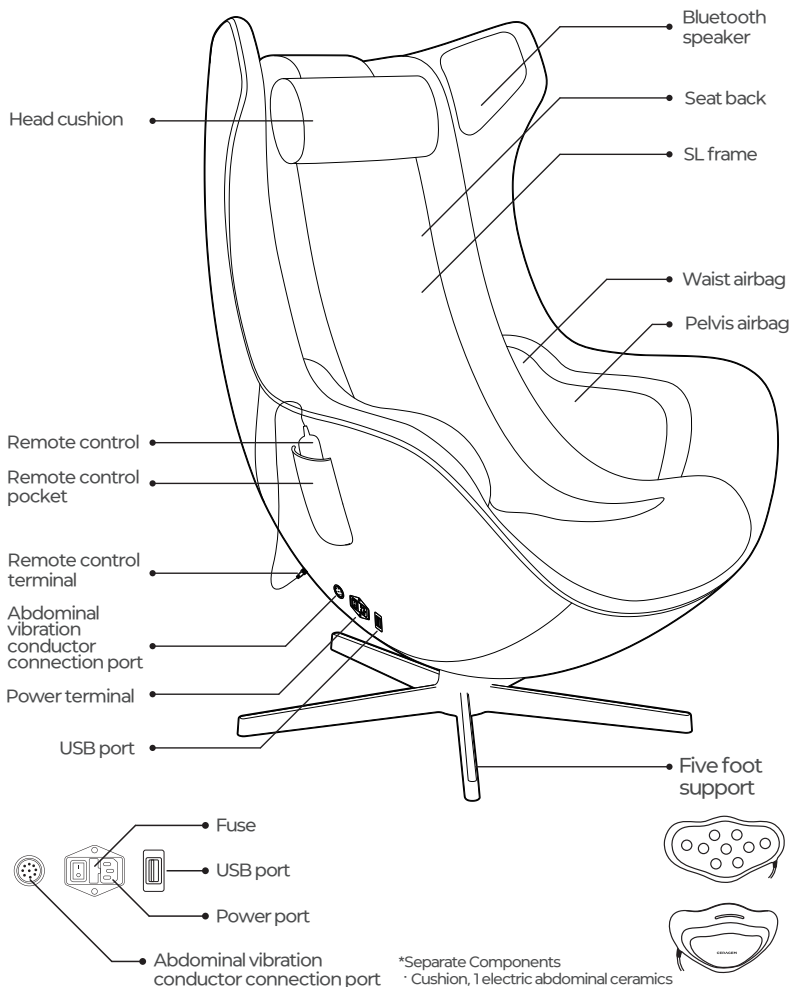

## 05 How to use the remote controller

---

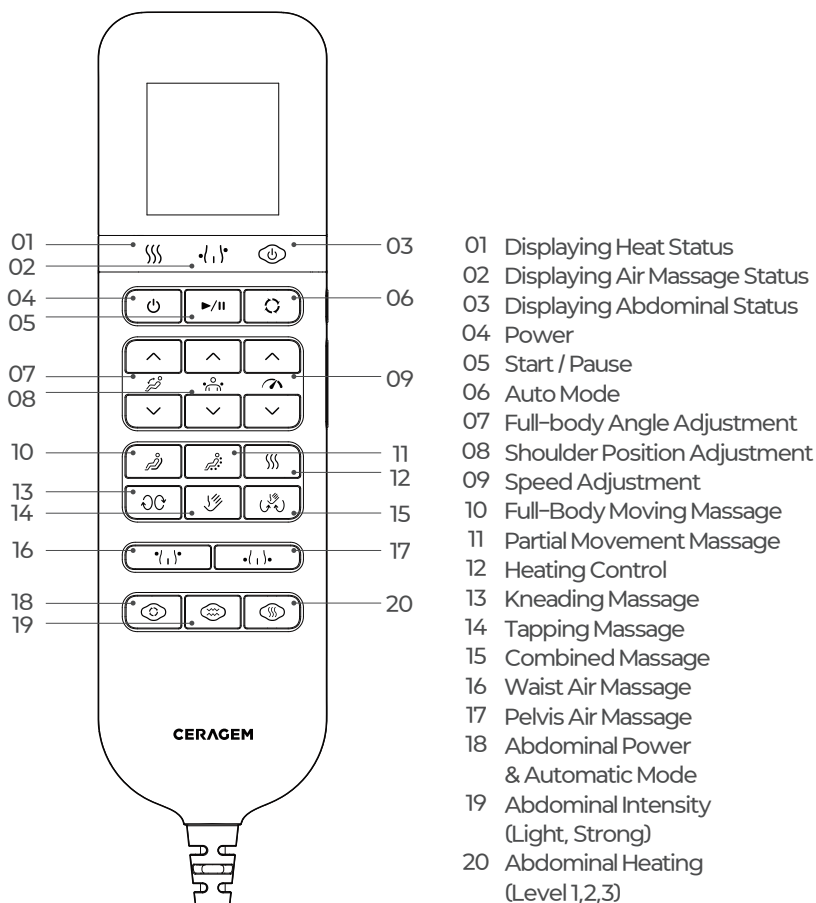



## Power

Power on and off the massage chair. If you press the 'power' button while in use, All massage movements stop, You can restore the original state by pressing the 'Start/Pause' button.

## Start / Pause

A button that allows you to 'start' or 'pause' the massage chair motion.

## Auto Mode

The buttons automatically implement the kneading/tapping/ combined massage in 8 modes. If you want to change to another mode while in use, you can set it up quickly by pressing Automatic Mode.

## Full-body Angle Adjustment

A button to adjust full body angle up/down.

## Shoulder Position Adjustment

The button to adjust the shoulder position in three stages: Up/Basic/Down.

## Speed Adjustment

The button to adjust the massage speed.

## Full-Body Movement Massage

Press the button to enable the massage ball to move over the entire body along the frame.

## Partial Movement Massage

Press the button to enable the massage ball to move within a certain range or a fixed range.

## Heating Control

A button that allows you to apply heat to the back and hips in three stages, lightly/medium/strongly, or turn off.

## Kneading Massage

It's a massage in which the massage balls gather to the side and then move away and rub them. The button can select forward/reverse/off by pressing sequentially.

## Tapping Massage

It is a massage button in which the upper and lower massage balls move back and forth. The button can select On/Off by pressing sequentially.

## Combined Massage

It is a massage button that repeats the complex massage method rhythmically. The button can select On/Off by pressing sequentially.

## Waist Air Massage

It's a button that allows you to get an air massage on your back. Can select light/strong/off by pressing sequentially.

## Pelvis Massage

It's a button that allows you to get an air massage on your Pelvis. Can select light/strong/off by pressing sequentially.

## Abdominal Power & Automatic Mode

This button that allows the automatic mode of the abdominal massager to be applied in mode 1 or mode 2. Can select mode1/mode2 /off by pressing sequentially.

## Abdominal Intensity

This button applies the vibration intensity of the abdominal massager in two stages: light, strong.

## Abdominal Heating

The button to select the thermal strength of the abdominal massager as light/medium/strong. You can select light/medium/strong/off by pressing sequentially.

## Auto Mode

Power-Auto Mode-Select the desired mode-Start.

### Auto Mode Start

1. ↻ Press Auto Mode to select the mode (eight).
2. ▶/⏸ Press the Start/Stop button to start reclining.
3. Once the reclining is complete, start the massage.








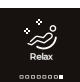
### Auto Mode Completion

1. At the end of the massage, the current reclining angle is maintained.
2. When the reclining notice appears, ▶/⏸ press the Start/Pause button to move to the original angle and exit.

### Power off during Auto Mode

1. Press the power button during the massage, the reclining guide will appear.
2. ▶/⏸ Press the start/pause button to back to original angle setting, shutting down the power.

### 8-Step Auto Mode

 Body concentration ■■■■■■■■■	Full-Body massage from neck to hip	 Neck & shoulders ■■■■■■■■■	Intensive massage from neck to shoulder
 Waist & Pelvis ■■■■■■■■■	Intensive massage from waist to pelvis	 Back & Waist ■■■■■■■■■	Intensive massage from back to waist
 Upper body Concentration ■■■■■■■■■	A massage that focuses from neck to waist	 Deep Tissue Massage ■■■■■■■■■	Full-Body complex massage with various roller directions
 Spine line ■■■■■■■■■	A combined massage that runs smoothly along the spine	 Relax ■■■■■■■■■	Relax massage to soothe tired body

※ If you select [Heat/Air Massage/Abdominal Massage] after the automatic mode massage starts, the function is added to the massage, and when the automatic mode ends, all are automatically shut down.

**Manual Mode**

Power – Choose from [full body/Partial Movement/Heating/Kneading/tapping /Combined Massage/air massage/abdominal massage] – Start.

**Massage Location****Manual Mode: Full Body/ Fix Range**

1.  Press the Manual – Full body button.
2. Press again at the desired position to fix the massage ball range.


**Manual Mode: Partial Movement/ Fix Range**

1.  Press the Partial movement/fix range button and select the Partial movement mode.
2. Timer remains when switching between other modes of operation (Kneading / Tapping / Combined) within manual mode operation.
3. Press one more time at the desired position to fix the massage ball.



**Power off during Manual Massage**

1. Press and hold the power button during the massage to get a reclining notice.
2.  Press the Start/Pause button to return to the original angle and off the power.

**Complete Manual Mode**

1. At the end of the massage, the current reclining angle is maintained.
2. When the reclining notice appears,  press the Start / Pause button to move to the original angle and exit.

**Manual Mode: Restart to another mode after completion**

1. At the end of the massage, the current reclining angle is maintained.
2. Press the Automatic mode or Manual mode button to select the desired mode.
3.  Press the Start/Pause button to start the massage.
4. If nothing is done at 1, it enters sleep mode.
5. Press any button to activate, then press the Auto or Manual mode button to select the desired mode.
6.  Press the Start/Pause button to start the massage.

※ If you select [Heating/Air Massage/Abdominal Massage] after the manual mode massage starts, the function is added to the massage, and when the manual mode ends, all are automatically shut down.

---

## Heat Massage

### Heating

1. Press button to adjust the temperature of the heat. [95°F/104°F/113°F]

※ When the heat function is turned on, the icon displaying the top status lights up.

## A Massage Method

### Kneading Massage

1. Press the kneading massage button to operate by selecting forward/reverse direction.
2. Press OFF to stop during the kneading massage.

### Tapping Massage

1. Press the tapping massage button to operate it.
2. Press OFF to stop during the tapping massage.

### Combined Massage

1. Press the complexing massage button to operate it.
2. Press OFF to stop during the complexing massage.

## Massage Speed

### Speed Adjustment

1. You can adjust the speed of the massage and the speed of movement up and down.
2. Press the speed control button to adjust the speed (low/medium/high).

## Additional Massage

### Air Massage: Waist

1. Press the button on the waist air massage to select strength (strong/light).
2. Press OFF to stop during the Waist air massage.

### Air Massage: Pelvis

1. Press the button on the Pelvis air massage to select strength (strong/light).
2. Press OFF to stop during the Pelvis air massage.

※ When the air massage is turned on, the icon displaying the top status lights up.

## Abdominal Vibration

### Power & Auto Mode

1. Ⓢ Press the power & Auto mode button of the abdominal pottery to operate by selecting the mode (1/2).
2. Press OFF to stop during the abdominal massage.

### Abdominal Intensity

1. Ⓢ Press the button on the abdominal massage intensity control button to select strength (light / strong).

### Abdominal Heating

1. Ⓢ Press the button on the abdominal heating control button to select strength (strong/normal/light).
  2. Press OFF to stop during the abdominal vibration.
- ※ When the abdominal massage function is turned on, the icon displaying the top status lights up.
- ※ If select multiple or consecutive other functions after the selected function starts, the function is added, and when the first selected function is terminated, all are automatically terminated.

## Bluetooth

### Bluetooth ON/OFF

1. Bluetooth ON: With OFF, press and hold the top button on the right side to turn on Bluetooth.
  2. Bluetooth OFF: With ON, press and hold the top button on the right side to turn off Bluetooth.
- ※ How to connect to Bluetooth: Turn on Bluetooth in your phone settings and locate and connect the device in the list.
- ※ Initial setting value.  
I D : CGMEMCG-2301  
PW : 1234

## Volume Control

1. Adjust the volume by pressing the + or - button located on the right side.
2. Press and hold the - button on the right side to mute.

If any problem occurs while using the product, refer to the following solutions.

If the problem persists, turn off the power supply, disconnect the power plug from the outlet, and then ask the service center.

Users should never disassemble or repair the product at their own discretion.

CERAGEM is not responsible for any damage to the product caused by disassembly or repair at the user's discretion.

Problem	Solution
The product does not work.	Check that the main power and power switch are turned on. Check the power plug connection state. Check that the correct voltage is used.
You can hear a fricative sound from the massage machine parts. You can hear the sound of the motor. Noise increases when using the tapping and kneading function.	Those noises are generated due to the mechanical structure and are not a malfunction. There may be noise even during normal operation.
The massage chair does not operate but the remote control does.	When the product is used for a long time, the massage chair automatically switches to overheating prevention mode due to increase in temperature. When the product is used for a long time, the massage chair automatically switches to overheating prevention mode due to increase in temperature. In this case, turn on the power and wait 30 to 50 minutes for the heat to cool.
Both the remote control and the massage chair do not work.	In case the working time has been set automatically, the set time has ended. Press the button again. Please turn the product off and on again.

**Q1 . Can I get a massage on my back only?**

After pressing the rubbing massage button, press the full-body moving massage button and press it again on the area where you want to be massaged. The massage ball will become fixed.

**Q2 . How far does the chair lean back?**

It reclines upto 145 degrees, from its default angle (115 degrees).  
The electric reclining function applied to the Pause M4 is adjustable.  
(It moves only when the user is pressing the angle adjustment button, and stops when not pressing.)

**Q3 I can't connect to Bluetooth Speaker.**

Turn the power off and then on. Then, delete CGMEMCG-2301 on the device to be connected, and reconnect it again.

\*Bluetooth can only be connected to one device.

**Q4 . Will it cause a lot of noise on the floor?**

Unlike other massage chairs, it barely causes noise on the floor because the main body is not attached to the floor.

In addition, you can place a rug under the massage chair to further reduce the noise.

**Q5 . Can I replace the leather?**

The head cushion, back seat, and air back seat can be replaced.  
The cost of replacement varies for each part.

**Q6 . The operating sound seems louder than what I heard when trying it at the store.**

At the store, you may not be able to hear the operating sound of the massage chair due to the ambient noise. In a relatively quiet residential space, you may hear the machine sound when operating the massage chair, but this is a normal and do not have to worry about it.

**Q7 . Is there a moving service for the chair when I move?**

Moving service is not provided and you must move the product by yourself.

**CERAGEM massage chairs are warranted for a period of ONE year against**

1. Missing materials from the Product at the time of delivery;
2. Defects in the manufacturing of materials used for the Product at the time of delivery;
3. Failure of the Product to correspond with the description in the applicable contract of sale at the time of delivery; and
4. Major defects in the functional department, such as :
  - Skeletal steel frame failure
  - Airbag failure
  - Power failure
  - Remote control failure
  - Massage ball failure
  - Recline tension mechanism failure

**Subject to the following**

Any replacement parts shall be sent to the address of the original purchaser, replaced parts may be required to be shipped back to CERAGEM at its discretion; to benefit under this warranty, CERAGEM must be notified of any defect within 14 days of the occurrence of such defect via an image or video of the relevant defect sent electronically to [ceragem.co.kr/customer/](mailto:ceragem.co.kr/customer/) and after inspection and acceptance by CERAGEM to the extent permitted under the applicable Federal and State law of the United States of America, all parts replaced or repaired under this warranty shall be warranted for the remainder of the designated Warranty Period.



### What Our Warranty Does Not Cover

All warranties listed on this page exclude the following where permitted under the applicable Federal and State law of the United States of America:

1. Normal wear and tear of surfaces that are not attributable to manufacturing, including but not limited to (a) degradation of upholstery (PU leather / fabric / genuine leather / mesh), top portion of armrests, pillows, casters exterior or mesh stitching, (b) damage to plastic coverings, (c) variance of tightness of mechanisms, and (d) excludes damage from accidents, natural disasters, physical damage from pets, fire and flood, and non-manufacturing incidents covered under insurance policies (regardless of whether you are in possession of the relevant insurance), (e) excludes any variations in dimensions and color;
2. Any claim for indirect, consequential, or incidental loss.
3. Damage to a product that has been modified or altered without the written permission of CERAGEM.
4. Damage caused to all CERAGEM products during self-assembly due to failure to adhere to the relevant self-assembly manual instructions provided;
5. Minor misalignments, variations arising out of hand-made workmanship. To ensure your product enjoyment and safety, and protection under the "Warranty Period", please ensure you have CERAGEM assembly and care guidelines for this product, and follow them for the correct assembly and care for your item. Contact our Customer Service Team to obtain these guidelines if you do not have them.

<b>Product Name</b>	MASSAGE CHAIR
<b>Model Name</b>	CGMEMCG-2301
<b>Size of the massage chair</b>	W29.5 X D44.1 X H39.4 in
<b>Size of the ottoman</b>	W20.1 X D16.1 X H18.5 in
<b>Net Weight</b>	83.1 lbs(massage chair),11 lbs (ottoman)
<b>Gross Weight (including box)</b>	130.9 lbs
<b>Working Voltage</b>	AC100-240V, 50/60Hz
<b>Power Consumption</b>	160W
<b>Manufacture</b>	Henglin Home Furnishings Co.,Ltd
<b>Country of Manufacturer</b>	China

**Wait!**

Warranty details can be found on [www.ceragemus.com](http://www.ceragemus.com)



FCC ID 2BBYO-CGMEMCG-2301

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

For warranty information,  
please visit the Warranty page at [www.ceragemus.com](http://www.ceragemus.com)

Product Information	
Product Name	Product Name
Manufacture Date . . .	Manufacture Number
Purchase Store	Purchase Store Contact
Date of Purchase . . .	Order Number
Date of Installation . . .	Warranty Date . . .
Customer Information	
Name	
Contact	
Address	

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This device complies with part 15 of the FCC Rules.

Operation is subject to the following two conditions:

1. This device may not cause harmful interference.
2. This device must accept any interference received, including interference that may cause undesired operation.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

#### FCC 20cm Statement

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This equipment should be installed and operated with a minimum distance of 20cm between the radiator and your body. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

NOTE

---

NOTE

---

# **CERAGEM®**

US Office  
CERAGEM International Inc.  
9350 Wilshire Blvd Unit 400,  
Beverly Hills, CA 90212

Customer  
Service Center  
9350 Wilshire Blvd Unit 400,  
Beverly Hills, CA 90212  
213 - 480 - 7070

# CERAGEM

US Office  
CERAGEM International Inc.  
9350 Wilshire Blvd Unit 400,  
Beverly Hills, CA 90212  
Copyright © All rights reserved by CERAGEM